

Rack Adapter Panel

GRA-446



Assembly Manual

GW INSTEK PART NO. 82RA-44600M01

This assembly manual contains proprietary information, which is protected by copyrights. All rights are reserved. No part of this manual may be photocopied, reproduced or translated to another language without prior written consent of Good Will company.

The information in this manual was correct at the time of printing. However, Good Will continues to improve products and reserves the rights to change specification, equipment, and maintenance procedures at any time without notice.

Good Will Instrument Co., Ltd.
No. 7-1, Jhongsing Rd., Tucheng Dist., New Taipei City 236, Taiwan.



ISO-9001 CERTIFIED MANUFACTURER

GW INSTEK

GRA-446


The GRA-446 rack mount is suited for the GPT-9500. The GRA-446 supports the EIA standards for use with standard-sized rack units. The GRA-446 holds one GPT unit and fits in a standard 19" 4U rack.

Rack Mount Adapter Component List and Assembly


Description	Please check to make sure all the components of the GRA-446 rack mount kit are provided.			
GRA-446 Assembly	Item	Part NO.	Description	QTY
	1	62RA-446HP101	L type bracket	1
	2	594B-W3008NS1	Screw BMS M3*8 ISO N (WASHER P.S), ROHS	3
	3	594B-W4008NS1	SCREW BMS ,+ ,4*08*0.7P ,ISO ,N (WASHER P.S) ,RoHS	5
	4	62RA-446LO101	GRA-446 Bottom cover	1
	5	62RA-446FP101	GRA-446 Front plate	1
	6	62PA-340HD101	CA PA 3U Handle, ROHS	2
	7	592B-W4016NS1	SCREW PMS ,+ ,4*16 ,ISO ,N , (WASHER P.S) ,RoHS	4

Assembly


- Initial setup
1. Remove a screw from (back front) the unit. 3 screws in total.

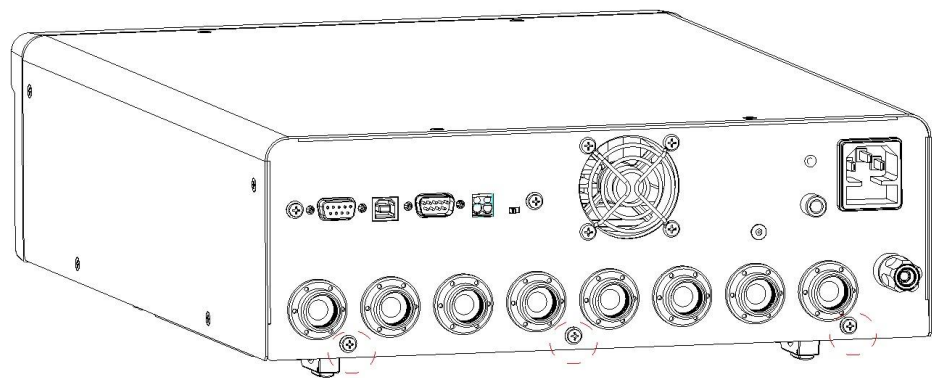

2. Attach the L-type brackets using the 594B-W3008NS1 screws to the unit.

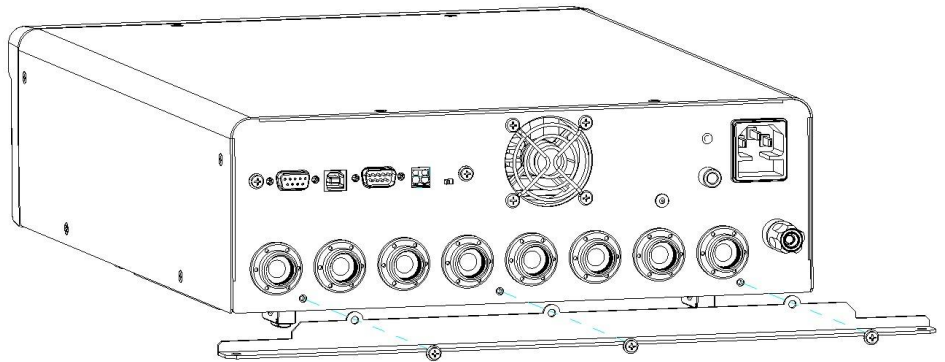

3. Attach the front plate and handles to the bottom cover using the screws.

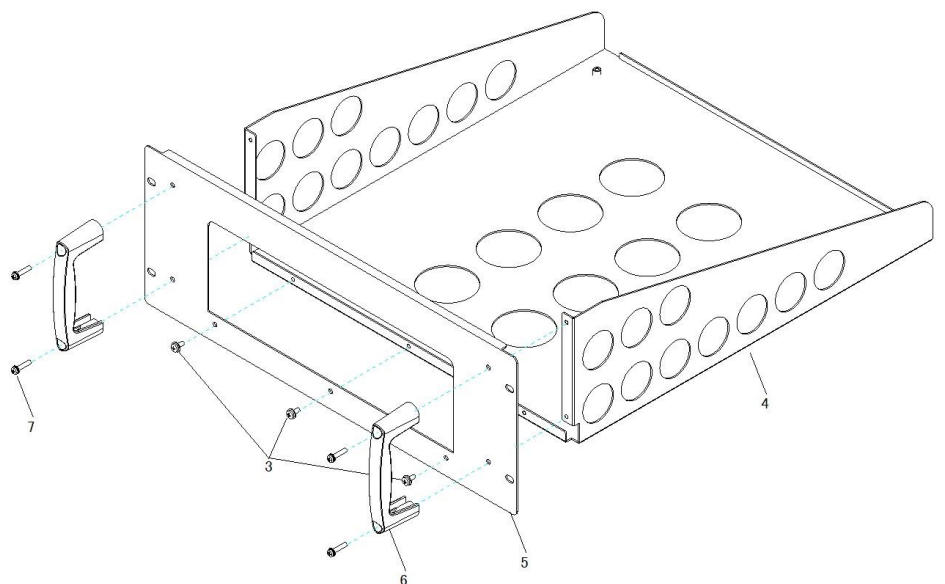

4. Attach the unit to the bottom plate with the L-brackets and the 594B-W4010NS1 screws. Make sure the front panel butts up against the front plate of the rack.

